

## Plate and Frame Filter BECO COMPACT™ PLATE 200



### Plate and Frame Filtration System for Laboratories and Research Centers

BECO COMPACT PLATE 200 depth filtration system from Eaton's Begerow Product Line is for small volume filtration tasks in laboratories and research centers. Applications include small batches and up-scaling in the pharmaceutical and biotechnology, chemical, cosmetics, and food industries.

BECO COMPACT PLATE 200 is a premium multi-sheet filter with a high compression of the depth filter sheets.


#### Special Features

- **Highly flexible** filter area (from 0.8 to 1.1 ft<sup>2</sup> (0.07 – 1.02 m<sup>2</sup>)) allows for customizable solutions for individual applications
- **Variable application options** for the use of filter plates, feed plates or cake frames. The filter can be used as a sheet filter, pre-coating filter or as a combination of both.
- **Easy handling** and **good cleanability** thanks to **ergonomic** design, smooth surfaces and open support grid of filter plates
- **Low drip loss** due to the compact design of the filter and a larger spindle
- **High chemical resistance** for the use of quality stainless steel materials, AISI 316L/AISI 316Ti

**Other designs** and further components for expansion to complete systems are available depending on individual requirements.

#### Technical Data

<b>Operating pressure:</b>	max. 145 psi (1,000 kPa/10 bar)
<b>Differential pressure:</b>	max. 58 psi (400 kPa/4 bar)
<b>Operating temperature:</b>	max. 284 °F (140 °C), depending on the filter plates used
<b>Filter area per sheet:</b>	0.37 ft <sup>2</sup> (0.034 m <sup>2</sup> )

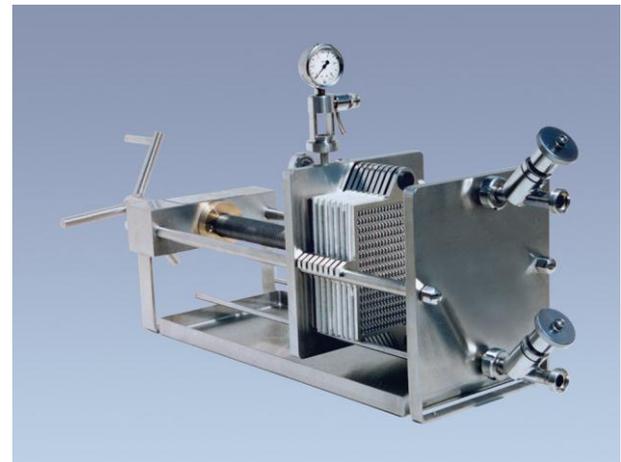

#### Materials

<b>Product contacting parts:</b>	Stainless steel AISI 316L/AISI 316Ti
<b>Other parts:</b>	Stainless steel AISI 304
<b>Gaskets/O-Rings:</b>	Silicone, Viton

#### Filter Plates/Cake Frames

	Width	Cake volume	Weight
Filter plate	0.32 in (8 mm)		1.8 lb (0.8 kg)
Feed plate	0.32 in (8 mm)		1.8 lb (0.8 kg)
Cake frame	0.32 in (8 mm)	0.03 gal (0.12 l)	1.8 lb (0.8 kg)
Cake frame	0.98 in (25 mm)	0.17 gal (0.65 l)	5.3 lb (2.4 kg)

#### Connections

<b>Inlet:</b>	DN 15, DIN 11851
<b>Outlet:</b>	DN 15, DIN 11851


#### Other connections

Flange according to DIN 2633/ISO or e.g. tri-clamp on request

Type	Filter area		Cake volume*		Filter plates/ cake frames	Depth filter sheets 7.9 x 7.9 in (200 x 200 mm)	Weight
	ft <sup>2</sup>	(m <sup>2</sup> )	gal	(l)			
200 SF	10.98	(1.02)	-		14/15	30	176 (80)
200 ASF with 0.32" (8 mm) cake frames	10.98	(1.02)	0.48	(1.8)	14/15	30	172 (78)
200 ASF with 0.98" (25 mm) cake frames	8.07	(0.75)	1.89	(7.15)	10/11	22	201 (91)

\* useable cake volume

Item	Description	Article no.	
1	Frame with tub	59262102	
2	First plate, 0.32" (8 mm)	59262301	
3	Filter plate, 0.32" (8 mm)	59262331	
4	Feed plate, 0.32" (8 mm)	59262332	
5	Cake frame, 0.32" (8 mm)	59262305	
6	Cake frame, 0.98" (25 mm)	59262304	
7	End plate, 0.32" (8 mm)	59262302	
8	Vent unit	59262312	
9	Inclined seat valve DN 15	54040204	
10	Gasket to items nos. 2, 3, 4, 5, 7	53126611 Silicone	53126613 Viton
11	Gasket to item no. 6	53126711 Silicone	53126713 Viton
12	Gasket DN 15, DIN 11851	50000152 Silicone	50000154 Viton

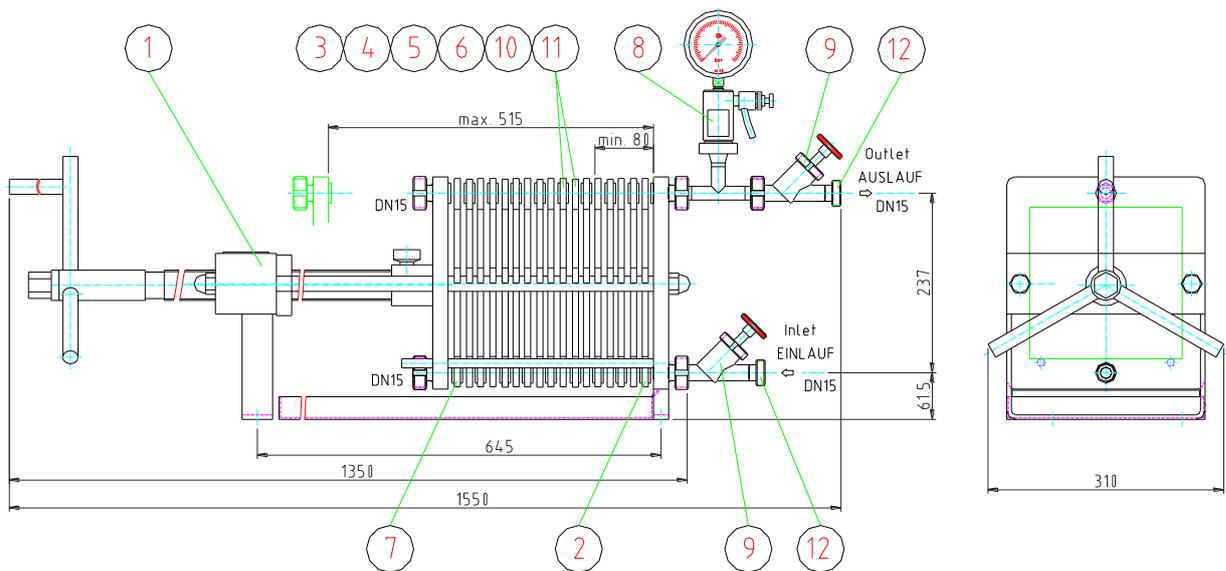



**North America - HQ**

44 Apple Street,  
Tinton Falls, NJ 07724  
Toll Free: 800 656-3344  
(North America Only)  
Voice: +1 732 212-4700

**Europe/Africa/Middle East**

Auf der Heide 2  
53947 Nettersheim, Germany  
Voice: +49 2486 809-0

**Internormen Product Line**

Friedensstraße 41  
68804 Altlussheim, Germany  
Voice: +49 6205 2094-0

**Begerow Product Line**

An den Nahewiesen 24  
55450 Langenlonsheim, Germany  
Voice: +49 6704 204-0

**Brazil**

Av. Julia Gaioli, 474 - Bonsucesso  
07251-500 - Guarulhos  
Brazil  
Voice: +55 11 2465 8822

**China**

No. 7 Lane 280 Linhong Road,  
Changning District,  
Shanghai 200335, China  
Voice: +86 21 5200 0422

**Singapore**

4 Loyang Lane #04-01/02  
Singapore 508914  
Voice: +65 6825 1668

For more information, please e-mail us at [filtration@eaton.com](mailto:filtration@eaton.com) or visit us online at [eaton.com/filtration](http://eaton.com/filtration) for a complete list of Eaton's filtration products.

Not all products in Eaton's Begerow Product Line are available in all regions. Please contact your local Eaton office to determine availability.

© 2012 Eaton Corporation. All Rights Reserved. All trademarks and registered trademarks are the property of their respective owners.

All information and recommendations appearing in this brochure concerning the use of products described herein are based on tests believed to be reliable. However, it is the user's responsibility to determine the suitability for his own use of such products. Since the actual use by others is beyond our control, no guarantee, expressed or implied, is made by Eaton as to the effects of such use or the results to be obtained. Eaton assumes no liability arising out of the use by others of such products. Nor is the information herein to be construed as absolutely complete, since additional information may be necessary or desirable when particular or exceptional conditions or circumstances exist or because of applicable laws or government regulations.



*Powering Business Worldwide*